

Technical Data Sheet

Doc No – Ab/TDS/GLP-110-II

Rev No : 02

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

TPX - Beakers without Handle

Valid for Catlog No. P50701, P50702, P50703, P50704, P50705, P50706, P50707, P50708

1. Description/Specification

1.1 Description:

Reusable TPX – Polymethylpentene (PMP) Beaker without Handle

1.2 Design:

Tapered Stacking Design

1.3 Performance Specifications:

1.3.1 Volume : 50 ml, 100 ml, 250 ml, 500 ml, 1000 ml, 2000 ml, 5000 ml, 10000 ml

1.3.2 Tolerance limit \pm 0.1 %

1.4 Material:

1.4.1 TPX- Polymethylpentene (PMP)

1.4.2 Comply with requirement of Directive EU No 10/2011 and requirement of FDA Regulation 21 CFR 177.1520 (c) 3.3b Polymers.

1.4.3 Materials are manufactured according to the requirements of GMP Regulation (EC) No 2023/2006

1.4.5 Material does not contain any Hazardous chemical substance & compliance with REACH AVHC & RoHS Free.

1.5 Color:

Translucent (Natural)

1.6 Graduations:

Large, Printed with temperature, chemical and solvent resistant ink.

1.7 Sterilization:

Can be autoclaved at 15 Psi, 121 °C for 20 minutes



Shiv Ganga Industrial Estate,
Plot no. 1 Lakeshwari, Bhagwanpur,
Roorkee - 247667 Uttarakhand, India
An ISO 9001:2008 Company
A Division of Abdos

Regd. Office :
13A, Government Place (East)
Kolkata – 700069, India
CIN : U85100WB2008PTC123903

P : +91 (0) 1332 235277
E : info@abdoslaware.com
W : www.abdoslaware.com

Technical Data Sheet

1.8 Quality control:

Testing of various characteristic in accordance with the valid Specification

1.9 other Information:

Reusable, can be autoclaved up to 10 autoclave cycle.

2. Basic Features

2.1 Made up of High purity Virgine Polymethylepentene (PMP)

2.1 Temperature range: -80°C to +175°C


2.3 Autoclavability: Yes, Autoclavable at 15 psi and 121°C for 20 minutes

2.4 Chemical resistance: See website.

Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121°C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50°C before removing product.


Quality Assurance Incharge
Abdos Labtech Pvt. Ltd.

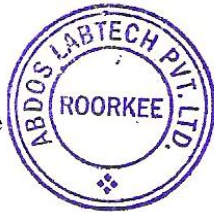

Shiv Ganga Industrial Estate,
Plot no. 1 Lakeshwari, Bhagwanpur,
Roorkee - 247667 Uttarakhand, India
An ISO 9001:2008 Company
A Division of Abdos

Regd. Office :
13A, Government Place (East)
Kolkata – 700069, India
CIN : U85100WB2008PTC123903

P : +91 (0) 1332 235277
E : info@abdoslabware.com
W : www.abdoslabware.com